

NANOpure® Diamond™ TOC Life Science (UV/UF with TOC) *ultrapure water system*

Operation Manual
Series 1197 and 1244

<u>Model No.</u>	<u>Voltage</u>
D11971	100V - 240V
D12441	100V - 240V

Table of Contents

Safety Information	3
Alert Signals.....	3
Warnings	3
Introduction	5
General Usage.....	5
General Specifications	6
Environmental Conditions	7
Declaration of Conformity	7
Unpacking	8
Installation	9
Choosing a Site	9
Ultrafilter Installation	10
UV Lamp Installation	11
Other Accessories	11
Bench Mounting	11
Wall Mounting	12
Installing the Control Panel in a Remote Location	12
Water Connections	13
Controls	15
Main Power Switch	15
Control Panel	15
Switches.....	16
Initial Operation	17
Cartridge Pack Installation.....	17
Initial Rinse	18
System Sanitization Procedure	19
Normal Operation.....	21
Dispensing Water.....	22
Automatic Dispensing	22
User Settings	23
Adjusting Display Brightness	23
Setting the Displayed Language	23
Setting the Date and Time	24
Setting the Cell Constant.....	25
Use of Standby Mode	26
Selecting the Set Point	26
Resetting the Sanitization Timer	27
Resetting the UV Timer	27
Resetting the TOC Lamp Timer	28
Using the TOC Analyzer	28
Setting Volumetric Dispensing	29
Setting Timed Dispensing	29
Performing a System Flush	30
Changing Purity Units	30
Temperature Compensation.....	30
Installing Float or Pressure Switch	31
N.I.S.T. Calibration Module	32
Performing an Electronic Calibration Using the Optional N.I.S.T. Calibration Module	32
Remote Dispenser	33
Computer/Printer Setup	34
Connecting NANOpure Diamond to Computer and Communicating Through the RS-232 Port Using Hyperterminal or Procomm	34
Connecting and Starting the Printer	36
Printer Setup	37
Maintenance and Servicing	39
General Cleaning Instructions	39
Cell Cleaning	39
System Sanitization	40
System Depressurization	41
Cartridge Pack Replacement.....	42
0.2 Micron Filter Replacement.....	43
Ultrafilter Replacement	43
UV Lamp Replacement	45
TOC Lamp Replacement	46
Fuse Replacement.....	47
Shutdown	47
Troubleshooting	48
Replacement Parts	51
Flow Chart Operational Overview	53
Wiring Diagram	54
Ordering Procedures	55
One Year Limited Warranty	56

Safety Information

Alert Signals



Warning

Warnings alert you to a possibility of personal injury.



Caution

Cautions alert you to a possibility of damage to the equipment.



Note

Notes alert you to pertinent facts and conditions.



Note

The UV lamp contains mercury. If broken or no longer needed, do not dispose of the UV lamp in the trash. Recycle or dispose of the UV lamp as hazardous waste.

Your Barnstead|Thermolyne NANOpure Diamond TOC Life Science (UV/UF with TOC) ultrapure water system has been designed with function, reliability, and safety in mind. It is your responsibility to install it in conformance with local electrical codes. This manual contains important safety information. You must carefully read and understand the contents of this manual prior to the use of this equipment. For safe operation, please pay attention to the alert signals throughout the manual.

Water purification technology employs one or more of the following: chemicals, electrical devices, mercury vapor lamps, steam and heated vessels. Care should be taken when installing, operating or servicing Barnstead products. The specific safety notes pertinent to this Barnstead product are listed below.

Warnings

To avoid electrical shock, always:

1. Use a properly grounded electrical outlet of correct voltage and current handling capacity.
2. Do not locate the NANOpure Diamond TOC Life Science (UV/UF with TOC) directly over equipment that requires electrical service. Routine maintenance of this unit may involve water spillage and subsequent electrical shock hazard if improperly located.
3. Replace fuses with those of the same type and rating.
4. Disconnect from the power supply prior to maintenance and servicing.

To avoid personal injury:

1. Do not use in the presence of flammable or combustible materials; fire or explosion may result. This device contains components which may ignite such materials.
2. This device is to be used with water feeds only. Sanitizing/cleaning agents must be used in compliance with instructions in this manual. Failure to comply with the above could result in explosion and personal injury.

SAFETY INFORMATION

3. Avoid splashing disinfecting or weak acidic solutions on clothing or skin.
4. Ensure all piping connections are tight to avoid chemical leakage.
5. Ensure adequate ventilation.
6. Carefully follow manufacturer's safety instructions on labels of chemical containers and material safety data sheets.
7. Depressurize system prior to removing the cartridge pack.
8. This unit is equipped with an ultraviolet lamp. Ultraviolet radiation is harmful to the eyes and skin. Do not attempt to observe the lamp directly.
9. Depressurize system prior to removing the ultra-filter.
10. Refer servicing to qualified personnel.

Introduction

Congratulations on your purchase of a Barnstead|Thermolyne NANOpure Diamond TOC Life Science (UV/UF with TOC) ultrapure water system. This water purification system is designed to provide nuclease-free (RNase and DNase), DNA-free, pyrogen-free, low TOC, high resistivity, reagent grade water that exceeds ASTM Type I, CAP and NCCLS Type I standards. It uses a four-stage deionization process combined with a UV lamp, an ultrafilter and a 0.2 micron filter to polish suitable feed water (distilled, deionized, or reverse osmosis) to produce pyrogen-free (<0.005 Eu/ml), low TOC (<3.0 ppb) water with a resistivity of up to 18.2 megohm-cm. Water resistivity is continuously monitored by a resistivity cell and displayed on a digital display. Actual system TOC values, obtained by an integrated TOC analyzer, can be concurrently displayed with the resistivity.

The NANOpure Diamond TOC Life Science (UV/UF with TOC) is designed to automatically clean the ultrafilter by flushing water through the filter for one minute every 12 hours. The water is sent to drain.

The resistivity electronics can be verified and calibrated utilizing a N.I.S.T. Traceable Calibration module. See accessory ordering information.

Please read the instructions carefully to ensure that you receive maximum benefit from the NANOpure Diamond TOC Life Science (UV/UF with TOC). Also, please fill out and return the enclosed warranty registration card as it will help us assure you of proper warranty coverage.

General Usage

Do not use this product for anything other than its intended usage.

General Specifications

Dimensions and Clearance Requirements

Dimensions

13.5" W x 19.5" H x 17.0" D (34.3 x 49.5 x 43.2 cm)

Clearances

Sides - 9" (22.9 cm) minimum for servicing.

Above - 3" (7.6 cm) minimum for removal of the top cover.

Front - 4.75" (12.1 cm) minimum for opening the front door.

Feed Water Requirements

Types	RO, DI, distilled.
TOC	< 1.0 ppm.
Turbidity	1.0 N.T.U. maximum.
Pressure Range	Gravity feed to 100 psig (7kg/cm ²) maximum.
Temperature Range	4°C - 40°C (40-104°F)
TDS (CaCO ₃)	< 70 ppm
Silica	< 1 ppm
Silt	< 5% SDI

Product Water

Water Quality

Resistivity	> 18.0 MΩ-cm
TOC	< 3.0 ppb
Endotoxin	< 0.005 Eu/ml

Flow Rate

1.5 lpm maximum at minimum inlet feed water pressure 10 PSIG with a new final filter.

Volumetric Dispense

Accuracy: ± 5%

Repeatability: ± 3%

TOC Analyzer

Range	1 to 250 ppb
Accuracy	±1 ppb or ±15% (whichever is greater)
Cycle Time	≤ 3.5 minutes
Water Consumption	≤ 15 ml per cycle
Min. Water Resistivity	> 10 MΩ for specified accuracy

Electrical Requirements

The NANOpure Diamond TOC Life Science (UV/UF with TOC) is equipped with 2 power cords to be plugged into an electrical outlet of the appropriate voltage.

U.K. customers use cord, CRX100 and fuses for 240V installation.

Model D11971 100-240 VAC, 100 watts, 47-63 Hz, 1 phase

Environmental Conditions

Operating: 4°C - 40°C; 20% to 80% relative humidity, non-condensing. Installation Category II (over-voltage) in accordance with IEC 664. Pollution Degree 2 in accordance with IEC 664.

Altitude limit: 3,500 meters.

Storage: -25°C to 65°C; 10% to 85% relative humidity.

Declaration of Conformity

Barnstead|Thermolyne hereby declares under its sole responsibility that this product conforms with the technical requirements of the following standards:

EMC: EN 61000-3-2 Limits for Harmonic Current Emissions
EN 61000-3-3 Limits for Voltage Fluctuations and Flicker
EN 61326-1 Electrical Equipment for Measurement, Control, and Laboratory Use; Part 1: General Requirements

Safety: EN 61010-1 Safety Requirements for Electrical Equipment for Measurement, Control and Laboratory Use; Part I: General Requirements

per the provisions of the Low Voltage Directive 73/23/EEC, as amended by 93/68/EEC.

The authorized representative located within the European Community is:

Electrothermal Engineering, Ltd.
419 Sutton Road
Southend On Sea
Essex SS2 5PH
United Kingdom

Copies of the Declaration of Conformity are available upon request.

Unpacking

1. Remove the unit from its shipping container. Remove all contents carefully. **Ensure that the UV lamp, ultrafilter, feed and drain tubing, sanitization cartridge, wall bracket, accessory parts bag and power cords (see list below) are removed from the packaging materials before discarding.** Put the NANOpure Diamond TOC Life Science (UV/UF with TOC) on a bench.

<u>Item</u>	<u>Part No.</u>
UV Lamp	LMX13
Sanitization Cartridge	CMX25
Wall Bracket (Unit)	BC1190X12
Wall Bracket (Remote Display)	BC1190X10
Ultrafilter	FL1192X1

Accessory Parts Bag Includes:

Feed Water Tube, 3/8" O.D.	TU1119X7
Drain Tube, 1/4" O.D.	TU1190X12
240 V Power Cord w/Fuses	CRX70 (2, FZX54)
120 V Power Cord w/Fuses	CRX72 (2, FZX47)
Display Cable	WHX20
Blank Display	DL1190X18
Plug Adapter	CEX42

Also required for installation:

* Customer Supplied Fasteners to wall mount unit and remote display.

Installation



Caution

Wall composition, condition and construction as well as fastener type must be considered when mounting this unit. The mounting surface and fasteners selected must be capable of supporting a minimum of 150 lbs.(68 kg). Inadequate support and/or fasteners may result in damage to mounting surface and/or equipment. If you are unsure of mounting surface composition, condition and construction or correct fasteners, consult your building maintenance group or contractor.



Warning

Do not locate the NANOpure Diamond TOC Life Science (UV/UF with TOC) directly over equipment that requires electrical service. Routine maintenance of this unit may involve water spillage and subsequent electrical shock hazard if improperly located.

Do not use in the presence of flammable or combustible materials; fire or explosion may result. This device contains components which may ignite such materials.



NOTE *****

The outlet of a gravity feed storage reservoir must be above or at the same level as the inlet of the NANOpure Diamond TOC Life Science (UV/UF with TOC).

Choosing a Site

The NANOpure Diamond TOC Life Science (UV/UF with TOC) features a removable control panel display which allows the system to be mounted almost anywhere within 10 ft. (3.1m) of the unit. Use the wall bracket for wall mounted systems as a template to drill mounting holes. (The NANOpure Diamond TOC Life Science (UV/UF with TOC) does not include screws and fasteners for mounting.) Please refer to the “General Specifications” section for clearance requirements.

INSTALLATION



Note

Do not install the ultrafilter until the NANOpure Diamond TOC Life Science (UV/UF with TOC) is ready for use. If the unit sits for longer than 72 hours without use, it will become contaminated and may affect product water quality.



Note

The ultrafilter must be installed before the 0.2 micron filter. Please refer to "0.2 Micron Filter Replacement" under "Maintenance and Servicing" for instructions on how to replace the 0.2 micron filter.



Top side port with cap removed

Bottom side port with cap removed

Ultrafilter with fittings



Top U-shaped
Tubing
Connection

Velcro strap

Bottom connector

Ultrafilter Installation


® Velcro is a registered Trademark of Velcro Industries B.V.

Ultrafilter Installation

Remove the ultrafilter from its protective bag. Install the ultrafilter as follows:

1. Open the front (right side) door and right side door.
2. Remove restraint (twist tie) from the two tubes.
3. Locate the metal bracket and separate the Velcro® strap.
4. Remove the cap on the top side port of the ultrafilter (make sure the bottom cap on the bottom side port is secure).
5. Align the ultrafilter with the metal bracket (Side ports should face the back of the unit) and push the bottom fitting tightly into the bottom connector.
6. Wet the end of the U-shaped tubing with deionized water before insertion.
7. Insert side ports into the holes in the metal bracket.
8. Push the U-shaped tubing into the top ultrafilter connector, making sure the connection is tight.
9. Secure the ultrafilter with the velcro strap.
10. Connect tubing with luer fitting to the top side port.
11. Close and latch right side door and front (right side) door.

Top side port
tubing with
luer fitting

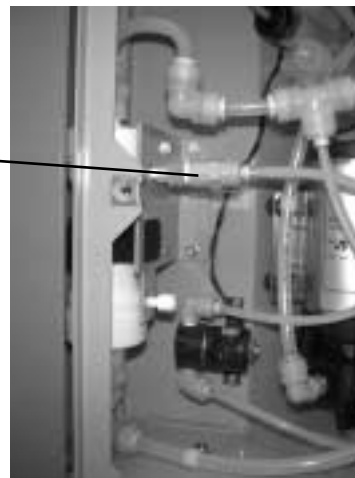

Ultrafilter Connections (Right side view)



Note

The UV lamp contains mercury. If broken or no longer needed, do not dispose of the UV lamp in the trash. Recycle or dispose of the UV lamp as hazardous waste.



Caution

DO NOT TOUCH THE GLASS PORTION OF THE LAMP! It is recommended that lint-free gloves be worn when handling the lamp. The glass portion must be free of fingerprints, perspiration, etc. Even a light coating of perspiration will reduce the effectiveness of the lamp. If the glass portion of the lamp is touched, clean it with a damp, lint-free cloth: use isopropyl alcohol as required.



UV Lamp Installation - Left side view

UV Lamp Installation

Locate the UV chamber inside the left door of the NANOpure Diamond TOC Life Science (UV/UF with TOC). Install the UV Lamp as follows:

1. Disconnect the power cord from the unit.
2. Access the left side of the unit by unlatching the screw securing the door with a flathead screwdriver. Open the door. The UV lamp will be installed vertically.
3. Remove black cap by sliding off UV housing. Slide the cap over the wires to allow the connector to be grasped.
4. Remove the UV lamp from its packaging. DO NOT TOUCH THE GLASS PORTION OF THE LAMP! It is recommended that lint-free gloves be worn when handling the lamp. The glass portion must be free of fingerprints, perspiration, etc. Even a light coating of perspiration will reduce the effectiveness of the lamp. If the glass portion of the lamp is touched, clean it with a damp, lint-free cloth: use isopropyl alcohol as required.
5. Insert the lamp halfway into the UV lamp chamber and plug it into the connector. Reposition the black cap and then fully insert the lamp into the chamber.
6. Make sure the UV lamp cable is tucked behind the chamber before closing and relatching the door.

Other Accessories

Optional accessories for this unit include a float switch, low pressure switch, N.I.S.T. Traceable Calibration Module, remote dispenser, and printer. Installation procedures for these accessories are explained in later sections of this manual.

Bench Mounting

1. Place the NANOpure Diamond TOC Life Science (UV/UF with TOC) on a bench top that is accessible to pretreated water, electricity and an atmospherically vented drain.

**Caution**

Wall composition, condition and construction, as well as fastener type, must be considered when mounting this unit. The mounting surface and fasteners selected must be capable of supporting a minimum of 150 lbs.(68 kg). Inadequate support and/or fasteners may result in damage to mounting surface and/or equipment. If you are unsure of mounting surface composition, condition and construction or correct fasteners, consult your building maintenance group or contractor.

**Note**

Please refer to the "General Specifications" section for clearance requirements.

**Warning**

Do not locate the NANOpure Diamond TOC Life Science (UV/UF with TOC) directly over equipment that requires electrical service. Routine maintenance of this unit may involve water spillage and subsequent electrical shock hazard if improperly located.

**NOTE *******

The outlet of a gravity feed storage reservoir must be above or at the same level as the inlet of the NANOpure Diamond TOC Life Science (UV/UF with TOC).

Wall Mounting

Install the NANOpure Diamond TOC Life Science (UV/UF with TOC) on a wall in a convenient location that is accessible to water, an atmospherically vented drain and electricity.

1. Locate the wall bracket packed separately from the unit.
2. Using the wall bracket as a template, locate and drill the mounting holes in the wall. A minimum of four (customer-supplied) fasteners will be required — two on the top and two on the bottom.
3. Attach the wall bracket to the wall using the customer-supplied fasteners.
4. Remove the locking screws on each side of the wall bracket.
5. Pull the two locking slides on each side of the wall bracket out as far as they will go.
6. Hang the unit on the wall bracket by sliding the mounting pins into the wall bracket slots.
7. Push the locking slides on each side of the wall bracket in as far as they will go.
8. Replace the locking screws.

Installing the Control Panel in a Remote Location

For your convenience, the control panel can be removed from the unit and mounted at a convenient location within 10 ft. (3.1 m.) of the unit. To remove the control panel from the unit:

1. Turn the unit OFF and disconnect it from the power supply.
2. Slide the control panel upward about 1/4" (.64 cm) to 3/8" (.95 cm) - (just enough to clear connector on top cover), pull forward, and remove. (The control panel fits tightly in the unit so it may be difficult to remove.)



Note

The removable control panel was not designed to be repeatedly removed from the unit, therefore, you may experience difficulty when attempting to remove it. It is therefore recommended that the control panel be permanently mounted in a remote location or remain intact with the unit.



Caution


Repeated removal and replacement of the control panel from the top cover may eventually cause it to become damaged.



Note

The blank protective cover on the rear side panel display connector should remain in place when the display is not remotely located.

Blank panel

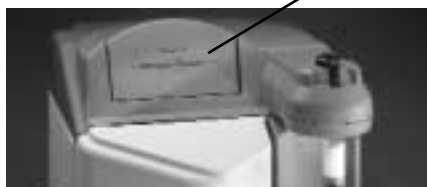

Control Panel
Wall Mount
Bracket



3. Replace the control panel with the blank panel provided with the unit to protect the electrical connector.
4. Using the wall bracket as a template, locate and drill the mounting holes in the wall. A minimum of two (customer-supplied) fasteners will be required. Attach the wall bracket to the wall using the customer-supplied fasteners.
5. Locate the 25-pin connector on the back top right of the unit. Remove the protective black cover and store it for future use if the control panel is returned to the top cover.
6. Locate the 10 ft. (3.1 m.) 25 pin M-F cable included with the unit; attach the female end to the 25 pin connector on the back top right of the unit, and the other end to the connector on the bottom of the control panel.
7. Tighten the cable screws on each end to secure the cable to the unit and control panel.
8. The control panel may now be mounted on the wall.

Water Connections

Feed Water Connection

1. Locate the length of 3/8" O.D. tubing provided with the unit with a quick disconnect insert on one end and a 3/8" (.95 cm) O.D. X 1/4" (.64 cm) NPT tubing adapter on the other.
2. Install the tubing adapter onto your incoming water line. We recommend a customer supplied shut off valve be installed in your feed water line. Do not connect the feed water to your NANOpure Diamond TOC Life Science (UV/UF with TOC). You will connect the feedwater during the "Initial Sanitization" procedure explained later in this manual. Note: The cartridge pack and ultrafilter must be installed before connecting water line.

INSTALLATION



Caution


Do not connect feed water until directed to do so during “Initial Sanitization.”

3. If a pressure switch is to be used, see the “Installing a Float or Pressure Switch” section for more information.

Atmospheric Drain Connection

When the NANOpure Diamond TOC Life Science (UV/UF with TOC) flushes its membrane, the water used is sent to drain through this connection. To install:

1. Locate the drain water tubing. This is the 1/4" (.64 cm) O.D. tubing that is approximately 6 ft. (1.9 m.) long. The atmospheric drain fitting is located on the lower left side of the NANOpure Diamond TOC Life Science (UV/UF with TOC).
2. Route the other end of the drain water tubing to an atmospherically vented drain and make a connection. Ensure there are no kinks in the tubing and that it proceeds in a downward plane. Proceed to “Initial Operation.”

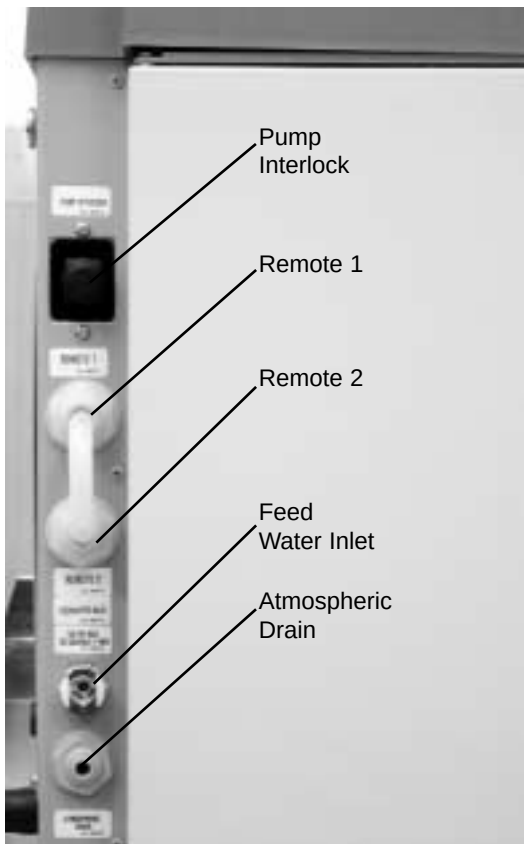

Water Connections - Left side view

Controls



Note

The removable control panel was not designed to be repeatedly removed from the unit, therefore, you may experience difficulty when attempting to remove it. It is therefore recommended that the control panel be permanently mounted in a remote location or remain intact with the unit.

Main Power Switch


The main power switch on the NANOpure Diamond TOC Life Science (UV/UF with TOC) is located on the lower back right side of the unit (as you face the front of the unit), directly above the power cord receptacle.

Control Panel

The NANOpure Diamond TOC Life Science (UV/UF with TOC) is controlled through a panel which incorporates switches to control its functions: START/STOP, STANDBY, DISPENSE, DOWN, UP, ENTER and BACK. This control panel utilizes a digital display to show system information and the purity of the product water. The intensity (brightness) of the display may be adjusted according to user preference. See “Adjusting Display Brightness” in the “User Settings” section of this manual.

A single LED located on the top of the control panel will illuminate to inform you of the system status as follows:

- LED is not illuminated: System is OFF or in (Idle) mode.
- Solid Green LED: System is recirculating and operating normally.
- Blinking Green LED: System is in Standby mode.
- Red LED: System is alerting user to an operational error or maintenance issue.

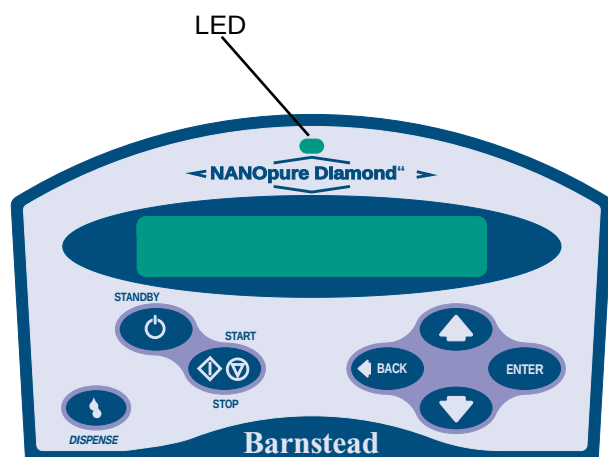

*NANOpure Diamond TOC Life Science
(UV/UF with TOC) Control Panel*

The control panel can be removed from the unit and remotely mounted. Please refer to “Installing the Control Panel in a Remote Location” in the “Installation” section of this manual for more information about how to remove the control panel from the unit.

Switches

When the main power switch (on the lower back right side of the unit as you face the unit) is ON, the switches on the control panel function as follows:

START/STOP allows you to alternate the unit between the normal (recirculation) mode; when the display is showing purity, and the (Idle) mode.

STANDBY allows you to put the unit into standby, recirculating water for 10 minutes/hour. The display will read “Standby” during periods of inactivity and “Recirculating” and the time remaining during the 10 minute recirculation.

DISPENSE allows you automatically deliver water from the unit. Please refer to “Automatic Dispensing” in the “User Settings” section of this manual for information.

UP and DOWN arrows allow you to scroll between menu options/items and values.

BACK allows you to return to the previous menu option/item. Please note that use of the BACK switch from a selection screen will return you to the previous menu option/item while maintaining the selection’s value upon entry.

ENTER allows you to activate a selected menu option/item and also functions as “Yes” whenever an option with a question mark appears.

DISPENSE KNOB when the unit is in the recirculation mode, push to the right to deliver water continuously until pushed back to the middle. Push to the left to deliver water manually until button is released.



Note

The unit **MUST** be in the normal recirculating mode (when the display is showing purity) in order to dispense water.

Initial Operation



Warning

Depressurize system prior to removing the cartridge pack.



Caution

Do not allow the NANOpure Diamond TOC Life Science (UV/UF with TOC) to operate unless water is available to the unit.



Note

Do not install the 0.2 micron filter and bell assembly at this time.



Cartridge Pack Types

The cartridge pack will come bagged with four manifold red protection caps.

Part Application

Diamond Kit Organic Free R/O & Distilled Feed

D50280 Ultra-Low Organics, Type 1 Water, Reverse Osmosis or Distilled Water Feed

Diamond Kit Organic Free Deionized Feed

D50281 Ultra-Low Organics, Type 1 Water, Deionized Water Feed

Each cartridge pack includes one 0.2 micron absolute final filter.

Cartridge Pack Installation

1. Open front (left side) door. Turn the cartridge pack so the red caps are facing you. Remove the red caps.
2. Verify that each of the four posts on the cartridge pack has an o-ring.
3. Lift unit manifold upwards, insert cartridge pack and align the pack so that it mates with the unit manifold.
4. Lower unit manifold until it is flush with the top of the cartridge pack.
5. Hand tighten wing head screws securely.
6. Close door.



Warning

Use a properly grounded electrical outlet of correct voltage and current handling capacity.

Initial Rinse

1. Attach the feed water line to the unit by snapping the quick disconnect coupling into the quick disconnect body in the lower left back of the unit.
2. Attach an atmospheric drain line [1/4" (.64 cm)] tubing by pushing into the quick connect fitting on the lower left back of unit. See layout connection under the "Water Connections" section.
3. Check to ensure the dispense knob is in the OFF (middle) position.
4. Select the appropriate power cord, remove and install the two fuses into the fuse drawer in the power entry module. U.K. customers use cord CRX100 and fuses for 240V installation.
5. Connect/attach to proper electrical outlet and the unit.
6. Turn the unit on by turning the main power switch to the "I" (ON) position.
7. The system greeting will display the type of unit you own. This greeting is "UV/UF with TOC."
8. The system electronics will initialize and check its calibration. (All units have been factory calibrated.) If calibration is OK the display will show "Calibration (passed)." The display will next read "Meter Self Test (in progress) Self Test. (passed)." Next, the display will read "10 (± 0.2) M Ω -cm." This is a reading of the calibration reference value. Finally, the display will show "TOC Initialization," "TOC: (In Progress)." (Note: The TOC initialization phase may take up to 45 seconds.)
9. From the (Idle) mode, when the display reads, "NANOpure Diamond (Idle) xx:xx:xx xx/xx/xx," press ENTER to run an air purge.
10. Unit will display, "Air Purge?" Press ENTER.
11. Upon completion of the air purge, press the UP or DOWN arrow until "TOC Short Flush?" displays. Press ENTER. (This flush will last approximately 20 seconds).

**Warning**

Avoid splashing disinfecting solution on clothing or skin.

Ensure all piping connections are tight to avoid leakage.

Ensure adequate ventilation.

Carefully follow manufacturer's safety instructions and material safety data sheets.

This device to be used with water feeds only.

Sanitizing/cleaning agents must be used in compliance with instructions in this manual. Failure to comply with the above could result in explosion and personal injury.

**Note**

After sanitization is complete, the sanitization timer will automatically be reset for approximately six months.



12. When the display returns to "TOC Short Flush?" press ENTER again to flush the TOC system a second time.
13. Press BACK to return to the (Idle) menu.
14. Complete the "System Sanitization Procedure."

System Sanitization Procedure

Your NANOpure Diamond TOC Life Science (UV/UF with TOC) has been shipped with a cleaning syringe (Cat. No. CMX25). When the system is new or whenever a new cartridge pack is installed, the system should be sanitized (cleaned). For UV/UF units, it is best to perform sanitizing procedures when the system or water is not required. The entire sanitization procedure will take approximately 3.5 hours to complete. Sanitize your NANOpure Diamond as follows:

1. From the (Idle) display, press the UP or DOWN arrow until display reads, "Sanitize Unit?"
2. Press ENTER.
3. Unit display will ask if you have "Disconnected Inlet Water (NO)."
4. Disconnect feed water supply at the quick-disconnect inlet valve on left side of unit.
5. Press the UP or DOWN arrow to select (YES) and then press ENTER. Display will now read, "Injected Sanitant (NO)."
6. Open front (left side door) to expose cartridge pack. Remove solution syringe from wrapper.
7. Remove luer cap on system injection port, located on the right side of the cartridge pack manifold by turning counterclockwise.
8. Remove syringe luer cap and attach syringe to the system luer fitting.
9. Slowly inject solution into system and remove syringe. Avoid injecting air into the system.
10. Replace luer cap on injection port and close door.

INITIAL OPERATION



0.2 Micron
Filter



Note

For more demanding applications where low TOC water is required, a rinse of 15-20 liters through the cartridges and 0.2 micron filter may be necessary.

11. Press the UP or DOWN arrow to select (YES) and then press ENTER. Display will now read, "Reconnected Water (NO)."
12. Attach feed water supply removed in step 4.
13. Press the UP or DOWN arrow to select (YES) and press ENTER. Display will now read, "Unit Sanitizing."
14. Unit may now be left unattended until display reads, "Sanitization Complete: Press ENTER."
15. After sanitization is complete the unit will return to the (Idle) mode.
16. Press the START/STOP switch to return the unit to the normal recirculation mode.
17. Remove the new 0.2 micron filter and bell assembly from the bag and insert it into the dispenser arm manifold. Gently turn it clockwise until it is fully seated.
18. Remove the protective cap from the filter bell. Turn the dispense knob to the left or right and flush 1-2 liters of water through the 0.2 micron filter. The product water is now ready for use. (See Normal Operation.)

Normal Operation



Note

On initial start-up, the purity meter may display "...MΩ-cm or uS/cm." This is caused by air in the cell and should change to a resistivity reading almost immediately. If the dots do not change to a value after one minute, refer to the "Troubleshooting" section of this manual.



Note

Do not turn off the NANOpure Diamond TOC Life Science UV/UF during non-work hours. Doing so will allow bacterial growth and other contamination of the water in the system. As a result, the system will require a rinse-up period at the beginning of the work day to achieve high quality product water. We recommend using the Standby Mode.



Note

From the normal recirculating mode when the display is showing purity, press ENTER to show the current temperature. The temperature will display for about 5 seconds before returning to the normal recirculating mode.



Note

Press any switch except ENTER to stop dispensing. Display will read, "Auto Dispense Interrupted by User."

1. Turn the system power ON by depressing the main power switch to the "I" position.
2. The system greeting display indicates the type of unit you own. This greeting is "UV/UF with TOC."
3. The system electronics will initialize and check its calibration. (All units have been factory calibrated.) If calibration is OK the display will show "Calibration (passed)." The display will next read "Meter Self Test (in progress) Self Test. (passed)." Next, the display will read "10 (±0.2) MΩ-cm." This is a reading of the calibration reference value. Finally, the display will show "TOC initialization," "TOC: (In progress)." (Note: The TOC initialization phase may take up to 45 seconds.)
4. From an (Idle) display, press the START/STOP switch on the control panel to enter the normal recirculating mode. The unit's pump will begin to run.
5. The display will then begin reading resistivity of the product water. Initially it will read "... MΩ-cm" until the air is bled from the system.
6. If the TOC analyzer has been chosen (See "Use TOC analyzer"), the display will begin showing TOC values of the product water. Initially it will read "TOC=...ppb" until valid data is available (which may take up to 4 minutes after starting).
7. Allow the water's resistivity to rise to the desired purity before drawing off water.
8. The system should be left on or in standby during the work day. (See "Use of Standby Mode.")

Dispensing Water

Use the dispense knob for manual dispensing or, use the “Dispense” switch for volumetric or timed dispensing. See “User Settings.”

Automatic Dispensing

(Factory set to “OFF.”) This feature is accessible only from the normal recirculating mode when the display is showing purity.

1. Place a container under the dispenser.
2. Press DISPENSE. **Water will start dispensing immediately!** Display will read, “xx.MΩ-cm xx Liters Remaining” or “xx.MΩ-cm xx Seconds Remaining,” or “AUTO DISPENSE <OFF>” depending on which dispensing method was last chosen in the auto dispense menu.

User Settings

Adjusting Display Brightness

1. From the (Idle) mode, when the display reads, "NANOpure Diamond (Idle) xx:xx:xx xx/xx/xx," press the UP or DOWN arrow until display reads, "Advanced Menu?"
2. Press ENTER.
3. Press ENTER to proceed through the options until "Set Display Brightness" is displayed.
4. Press the UP or DOWN arrow to choose "Yes."
5. Press ENTER.
6. Press the UP or DOWN arrow until the desired value (100%, 75%, 50% or 25%) is displayed.
7. Press ENTER.
8. Press BACK until you return to the (Idle) display.

Setting the Displayed Language

1. From the (Idle) mode, when display reads, "NANOpure Diamond (Idle) xx:xx:xx xx/xx/xx," press the BACK switch until the display reads "LANGUAGE."
2. The currently selected language is displayed on line 2 of the display.
3. Press the UP or DOWN arrow to select the language desired.
4. Press ENTER.

**Note**

Time values are displayed in 24 hr. military time. For example; 6:00p.m. = 18:00:00 = hr./min./sec.)

Setting the Date and Time

1. From the (Idle) mode, when the display reads, "NANOpure Diamond (Idle) xx:xx:xx xx/xx/xx," press the UP or DOWN arrow until display reads, "Advanced Menu?"
2. Press ENTER.
3. Press ENTER when the unit displays "Printer Attached?"
4. Continue pressing ENTER until "Date Format" is displayed.
5. When prompted with "Date Format" press the UP or DOWN arrow to show the date format desired. Press ENTER to select the desired format and continue.
6. When prompted with "Set Date," press the UP or DOWN arrow to choose "Yes," then press ENTER.
7. Press the UP or DOWN arrow until you reach the desired month. Press ENTER.
8. Press the UP or DOWN arrow until you reach the desired day. Press ENTER.
9. Press the UP or DOWN arrow until you reach the desired year. Press ENTER.
10. When prompted with "Set Time" press the UP or DOWN arrow to choose "Yes," then press ENTER. This is to set the time of day.
11. Press the UP or DOWN arrow until you reach the desired hour value. Press ENTER.
12. Press the UP or DOWN arrow until you reach the desired minute value. Press ENTER.
13. Press BACK twice to return to the (Idle) mode.

Setting the Cell Constant

The actual temperature and cell constants are attached to the cell cable internal to the NANOpure Diamond unit. This data helps ensure that the purity displayed is as accurate as possible. The only time the user would need to enter this data is if the cell or main board is being replaced and/or updated.

1. Shut power OFF to the unit and disconnect the power cord.
2. Unlatch and open the right side door.
3. Find and copy down the two cell constants attached to the gray cable, leaving the electronics board mounted on the inside of the right side door.
4. Close and latch the right side door.
5. Reconnect the power cord and turn unit ON.
6. From the (Idle) mode, when the display reads, "NANOpure Diamond (Idle) xx:xx:xx xx/xx/xx," press the UP or DOWN arrow until the display reads, "Advanced Menu."
7. Press ENTER.
8. Proceed through the options by pressing ENTER until "Set Cell Constant" is displayed.
9. Press the UP or DOWN arrow to select "Yes."
10. Press ENTER.
11. Note that the digits are entered from left to right.
12. Press the UP or DOWN arrow to change the "blinking" digit. Once the digit is correct, press ENTER to advance to the next digit.
13. After all digits have been entered, line 2 of the display will read, "(NO) Finished."
14. If the cell constants displayed on line 1 of the display are correct, press the UP or DOWN arrow to select "Yes."

15. Press ENTER.
16. Press BACK twice to return to the (Idle) mode.

**Note**

When in Standby mode, press START/STOP to continue normal operation.

Use of Standby Mode

At the end of the work day, place the NANOpure Dlamond TOC Life Science (UV/UF with TOC) in Standby mode for the night. To place the unit in Standby mode:

1. From the (Idle) mode, when display reads, "NANOpure Dlamond (Idle) xx:xx:xx xx/xx/xx," press STANDBY.

—OR—

2. From the xx.x MΩ-cm or x.xx uS/cm operational display, Press STANDBY.

After STANDBY is pressed the display will read "Standby" and the green LED will slowly blink. For 10 minutes of every hour while the pump is energized, the display will read "Standby (Recirculating) xx minutes." Every fourth time the pump goes on in the standby mode, the UV light will go on. The lamp is on 10 minutes out of every four hours. Press START/STOP to return to the recirculating mode.

Selecting the Set Point

The NANOpure Dlamond TOC Life Science (UV/UF with TOC) electronics include a user programmable set point which alerts you when water quality falls below the programmed set point. The set point is user selectable from 1-18 MΩ-cm (factory set at 10 MΩ-cm.).

1. From the (Idle) mode, when display reads, "NANOpure Dlamond (Idle) xx:xx:xx xx/xx/xx," press the UP or DOWN arrow until display reads, "Advanced Menu?"
2. Press ENTER.
3. Proceed through the options by pressing ENTER until "Adjust Set Point" is displayed.
4. Press the UP or DOWN arrow to choose "Yes."
5. Press ENTER.
6. Press the UP or DOWN arrow until the desired value is displayed.

7. Press ENTER.
8. Press BACK to return to the (Idle) mode.

The second line of the display will show a "Below Set point" message once every 8 to 10 seconds when the measurement is below the set point. (The set point indication is set at the factory at 10 MΩ-cm.)

Resetting the Sanitization Timer

1. From the (Idle) mode, when the display reads, "NANOpure Diamond (Idle) xx:xx:xx xx/xx/xx," press the UP or DOWN arrow until the display reads, "Reset Timers?"
2. Press ENTER. Display will read, "Reset UV Timer (NO)."
3. Press ENTER. Display will read, "Reset TOC Timer (No)."
4. Press ENTER.
5. Display will read, "Reset San Timer" and show on the right side of the second line when the user will start receiving "Sanitize Unit" reminders.
6. Press the UP or DOWN arrow to select "Yes."
7. Press ENTER and the timer will be reset. This will reset the sanitization timer for approximately six months calendar time.

Resetting the UV Timer

When the NANOpure Diamond TOC Life Science (UV/UF with TOC) is in an (Idle) position:

1. From the (Idle) mode, when display reads, "NANOpure Diamond (Idle) xx:xx:xx xx/xx/xx," press the UP or DOWN arrow until the display reads, "Reset Timers?" Press ENTER.
2. Press ENTER until the display reads, "Reset UV Timer." It will also read, on the right side of the second line, the hours of usage which have occurred since the UV Timer was last reset.

3. Press the UP or DOWN arrow to select "Yes."
4. Press ENTER and the timer will be reset for 1820 hours of operational time.

Press BACK two times to return to the (Idle) mode.

Resetting the TOC Lamp Timer

When the NANOpure Diamond TOC Life Science (UV/UF with TOC) is in an (Idle) position:

1. From the (Idle) mode, when display reads, "NANOpure Diamond (Idle) xx:xx:xx xx/xx/xx," press the UP or DOWN arrow until the display reads, "Reset Timers?" Press ENTER.
2. Press ENTER until the display reads, "Reset TOC Timer." It will also read, on the right side of the second line, the hours of usage which have occurred since the UV Timer was last reset.
3. Press the UP or DOWN arrow to select "Yes."
4. Press ENTER and the timer will be reset for 1300 hours of operational time.

The display will now display "Reset San Timer." Press BACK two times to return to the (Idle) mode.

Using the TOC Analyzer

(Factory set to "No.")

This feature is accessible from the (Idle) mode.

Use the analyzer when TOC values and resistivity/conductivity readings are to be displayed. DO NOT use the analyzer when only resistivity/conductivity readings are to be displayed.

1. From the (Idle) display press UP or DOWN arrow until "Advanced Menu?" is displayed.
2. Press ENTER.
3. Press ENTER when "Printer Attached?" is displayed.
4. Continue pressing ENTER until "Use TOC Analyzer?" is displayed.



Note

Volumetric dispensing is based upon total fluid volume entering the NANOpure Diamond TOC Life Science (UV/UF with TOC). Point of drawoff accuracy will be compromised if you are dispensing water volumetrically while simultaneously drawing water manually through the remote dispenser accessory.



Note

Once Volumetric or Timed Dispensing values are set, they will remain in memory until changed by the user. To turn the dispense method OFF, select "OFF" in the "Dispense Method" menu.

5. Press UP or DOWN arrow to choose "Yes."
6. Press ENTER.
7. Press BACK twice, from the "Date Format" display, to return to the (Idle) mode.

Setting Volumetric Dispensing

1. From the normal recirculating mode, when the display is showing purity, press the UP or DOWN arrow until the display reads, "Auto Dispense Menu?"
2. Press ENTER.
3. Display will read, "Dispense Method." Press the UP or DOWN arrow to select VOLUMETRIC.
4. Press ENTER. Display will read, "Volume Setting: xx Liters."
5. Press the UP or DOWN arrow to set the desired volume value. Press ENTER.

UP arrow: +0.25 up to 1 L, then 5 L, then incremental by 5 L

DOWN arrow: Decrement by 1 L until 1 L, then by 0.25 L

Maximum: 60 L

Minimum: 0.25 L

Setting Timed Dispensing

1. From the normal recirculating mode, when the display is showing purity, press the UP or DOWN arrow until the display reads, "Auto Dispense Menu?"
2. Press ENTER.
3. Display will read, "Dispense Method." Press the UP or DOWN arrow to select TIME.
4. Press ENTER. Display will read, "Time Setting: xx Minutes."
5. Press the UP or DOWN arrow to set the desired time value. Press ENTER.

UP arrow: +1 up to 5 min., then increments by 5 minutes
DOWN arrow: Decrement by 1 minute
Maximum: 40 minutes
Minimum: 1 minute

**Note**

An automatic 1-minute flush will occur once every 12 hours. If a unit is (Idle) or in STANDBY for over 12 hours, an automatic system flush will occur immediately after the START/STOP switch is pressed.

**Note**

Press the STOP/STOP switch to cancel the flush.

Performing a System Flush

1. From the normal recirculating mode, when the display is showing purity, press the UP or DOWN arrow until the display reads, "System Flush?" Press ENTER.
2. Press the UP or DOWN arrow to select "Short Flush?" or "Extended Flush?" Press ENTER.
3. Display will read, "System Flushing: (xxx) Seconds or (xxx) Minutes Remaining."

Changing Purity Units

(Factory set to "MΩ-cm.")

1. From the normal recirculating mode, when the display is showing purity, press the Up or DOWN arrow until display will reads, "Change Units?" Press ENTER.
2. Press the UP or DOWN arrow to select "MΩ-cm" or "uS/cm" Press ENTER.
3. Press BACK to return to the normal recirculating mode.

**Note**

Purity displayed will have an asterisk (*) the the left of the reading (* xx.x MΩ-cm or * x.xx uS/cm) to inform the user that uncompensated purity is being displayed.

Temperature Compensation


(Factory set to "ON.")

1. From the normal recirculating mode, when the display is showing purity, press the Up or DOWN arrow until the display reads, "Temperature Compensation?" Press ENTER.
2. Press the UP or DOWN arrow to select "OFF" or "ON" Press ENTER.
3. Press BACK to return to the normal recirculating mode.


Installing Float or Pressure Switch

Accessories D8964 (float switch) and D2706 (pressure switch) are designed to protect the NANOpure Diamond TOC Life Science (UV/UF with TOC) pump by alerting the system of an inadequate feed water condition so that the pump can be shut down. If an inadequate feed water condition exists and the pump is shut down due to this condition, the display will read "Check Inlet." Follow the proceeding instructions for installation.

1. Disconnect the unit from the electrical power.
2. If using the D8964 float switch, follow the installation instructions included with the float switch for installation into a tank.
3. If using the D2706 low pressure switch, install the PVC tee (supplied with D2706) in the incoming water line. Screw the switch into the top of the tee, then connect the inlet tubing of the NANOpure Diamond TOC Life Science (UV/UF with TOC) into the remaining opening of the PVC tee.
4. Route the cable from the float or low pressure switch to the top left of the NANOpure Diamond TOC Life Science (UV/UF with TOC).
5. Remove the jumper plug from the back, upper left portion of the unit and save for future use.
6. Plug the cable into the jumper plug outlet.
7. Reconnect the electrical power.

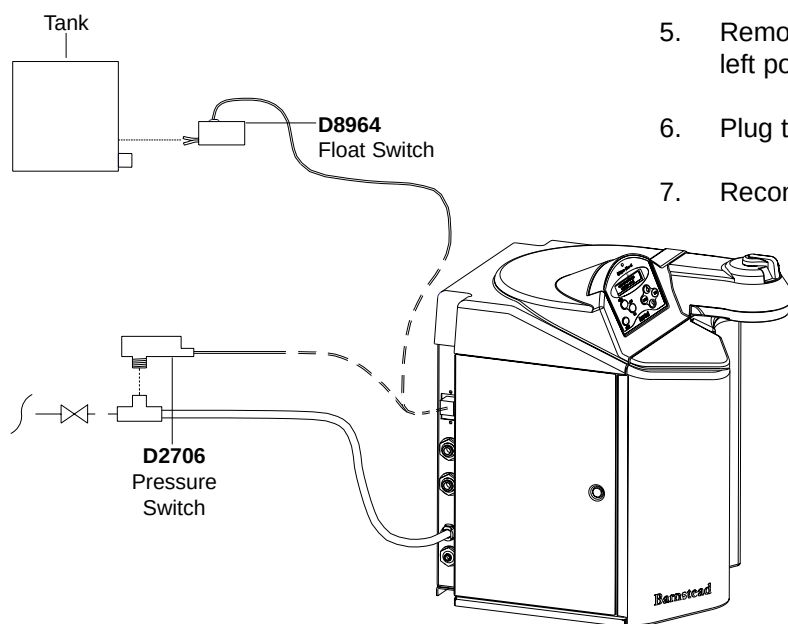

N.I.S.T Calibration Module

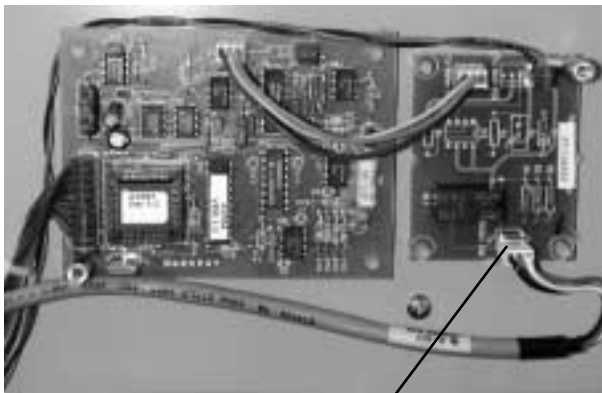

PC Board

J2 Connection



Note

The calibration (in progress) may take up to two minutes to complete.

Performing an Electronic Calibration Using the Optional N.I.S.T. Calibration Module

If you purchased the optional N.I.S.T. calibration module (Catalog No. E896X5) you can perform a calibration of the NANOpure Diamond purity sensing electronics traceable to N.I.S.T. standards.

1. Disconnect power cord from the unit.
2. Unlatch and open the right door. Remove the 2 screws and remove the top PC board cover plate over the PC boards.
3. Locate the outlet cell connection on the circuit board and remove the outlet cell connection at J2.
4. Connect the N.I.S.T. calibration module to the connection point (J2), "To Outlet Cell" on the calibration board (PC1190X2).
5. Leave the door open. Otherwise, the calibration module may collide into the UF plumbing.
6. Reconnect the power cord and turn the unit on by depressing the main power switch.
7. From the (Idle) mode, Press the UP or DOWN arrow until the display reads, "N.I.S.T. Calibrate?" Press ENTER.
8. Display will show "Calibration Module Installed?" Press the UP or DOWN arrow to select (YES).
9. Press ENTER.
10. Display will read "Calibration (In Progress)." When the calibration is complete the display will read "Calibration (passed/failed) for 4 seconds and then return to the "N.I.S.T. Calibrate?" display.

11. Press BACK to return to the (Idle) mode.
12. Turn the main power switch to the OFF "O" position. Unplug the unit.
13. Remove the calibration module and reconnect the cell.
14. Reattach the cover plate and close and latch the door.
15. Reconnect the power cord.

The calibration procedure is complete.

Remote Dispenser

Refer to the Operating Instructions for your remote dispenser.

Computer/Printer Setup

Connecting NANOpure Diamond to Computer and Communicating Though the RS-232 Port Using Hyperterminal or Procomm

Hyperterminal

RS-232 Capture Instructions

1. Connect the 9 pin serial cable (part no. WHX18) from the NANOpure Diamond RS-232 port to COM2: port or (COM1: port) on back of computer.
2. From the (Idle) display, press the UP or DOWN arrow until display reads, "Advanced Menu."
3. Press ENTER.
4. From the "Printer Attached" display, press the Up or DOWN arrow for (NO).
5. Press ENTER.
6. Now, during the normal Recirculation mode, date, time, resistivity/conductivity and temperature readings will be sent to the computer once per minute. If "Use TOC Analyzer" has been set to "Yes," the TOC data will also be sent.
7. Run your Microsoft Windows program. Open Hyperterminal (located in the accessories directory), give your setting a name and choose a symbol.
8. In the "Connect Using:" box, select "direct to" COM2: or (COM1:) as port (depending on your computer), and click OK.
9. Using the selection boxes in the next screen, configure the options as 9600 baud, 8-bit, No parity, 1 stop bit, Flow Control: None.
10. To receive data to a file, select Transfer and capture text. Accepting default will put the file



Note

Due to character set differences, hyperterminal may display the Ω symbol as $\hat{e}$. Also, the $^\circ$ symbol may be displayed as $\emptyset$.

which you name in the c:/windows directory called capture.txt and the file will be an ASCII text file.

11. You should see the data on the computer screen.
12. When you are finished storing data from the NANOpure Diamond, save your file. You can exit Hyperterminal or set up a new data file.
13. Hyperterminal data can be manipulated and graphed in Excel®.

Procomm

1. Connect the 9 pin serial cable from the NANOpure Diamond RS-232 port to COM2: port or (COM1: port) on back of computer.
2. From the (Idle) display, press the UP or DOWN arrow until the display reads, "Advanced Menu."
3. Press ENTER.
4. From the "Printer Attached" display. Press the Up or DOWN arrow for (NO).
5. Press ENTER.
6. The unit will now, during the normal Recirculation mode, send a resistivity/conductivity and temperature reading to the computer once per minute. If "Use TOC Analyzer" has been set to "Yes," the TOC data will also be sent.
7. Open DOS window and change directory to procomm directory.
8. Open procomm executable.
9. Press ALT-F10 for configuration screen.
10. Enter into Line Settings screen by typing ALT-P.
11. Type in 11 <Enter>, sets parameters as 9600, 8, N, 1.

12. Type in 21 <Enter>, sets to COM2: or (Type in 20 <Enter>, sets to COM1:), depending on your computer.
13. Type in 24 <Enter>, saves settings.
14. Type ESC to exit setup menu.
15. Procomm is now ready to accept input from the NANOpure Diamond.
16. To begin downloading a file, Type ALT-F1, and give the file an appropriate name when prompted.

The optional printer (part no. AY1137X1) provides a paper tape record for future reference.

During normal Recirculation, the printer prints date, time, purity and temperature readings and the optional TOC data, once every 12 minutes.

Connecting and Starting the Printer

1. Make sure the printer and NANOpure Diamond are turned OFF.
2. Connect the NANOpure Diamond to the printer via the RS232 port. Use the tan printer cable included with the Epson printer. Plug the printer power cord into an electrical outlet.
3. Turn the NANOpure Diamond ON.
4. Turn the printer ON.
5. From the (Idle) display, press the UP or DOWN arrow until display reads, "Advanced Menu."
6. Press ENTER.
7. From the "Printer Attached" display, press the Up or DOWN arrow for (YES).

8. Press ENTER.
9. Press BACK twice to return to the (Idle) mode.
10. During normal Recirculation mode, a resistivity/conductivity and temperature reading will be printed once every 12 minutes. If "Use TOC Analyzer" has been set to "Yes," the TOC data will also be printed.

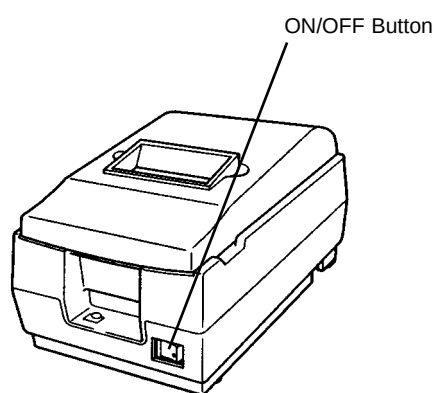


Figure A: Printer

Printer Setup

Power Switch

Located on the front of the printer (AY1137X1), this switch turns power to the printer ON and OFF. (See Figure A.) The green POWER light will illuminate when the printer is ON.

Paper Feed Button

Press this button once to advance the tape one line, or hold the button down to advance paper continuously. When the paper roll nears the end, the red PAPER OUT light will illuminate.

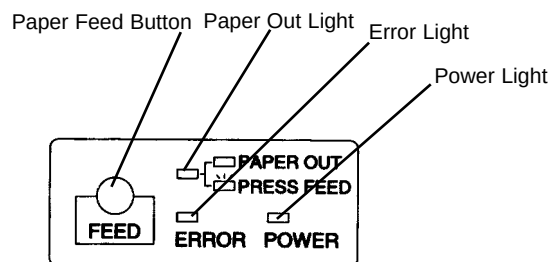


Figure B: Control Panel

If the red ERROR light illuminates, the printer is off line. This could be a result of the print head being too hot, in which case the printer will resume printing once it cools down; or—the paper may be jammed in the printer. Turn the printer OFF before checking for jammed paper. Remove the paper jam and turn the printer back ON. If the printer still will not print, unplug it and refer servicing to qualified personnel.



Note

A colored stripe along the edge of the paper indicates that you are approaching the end of the tape. When the stripe appears, approximately 23 inches (58.4 cm) of paper remains, enough for 1-2 days. Replace the tape before it is completely exhausted to ensure a complete record of your analysis.

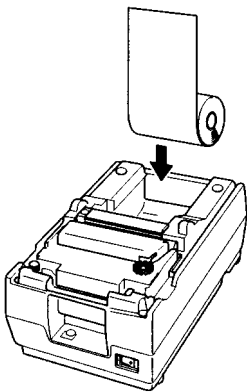


Figure C: Inserting the Paper Roll

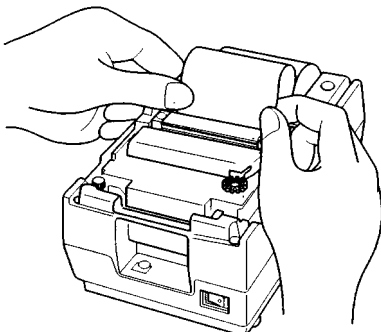


Figure D: Inserting the Paper Into the Paper Slot

Installing the Paper Roll

When the red PAPER OUT light illuminates, it is time to replace the paper roll. Replace the paper roll (AY669X4) as follows:

1. Turn the printer ON and open the chamber cover.
2. Make sure the edge of the paper roll is straight.
3. Insert the paper roll into the printer. (See Figure C.)
4. Insert the paper straight into the paper slot (see Figure D). The paper will feed automatically.
5. Tear off a strip of paper and close the cover. (The red PAPER OUT light will be OFF if paper roll installation is successful.)

Maintenance and Servicing



Warning

- Avoid splashing disinfecting or weak acidic solutions on clothing or skin.
- Ensure all piping connections are tight to avoid chemical leakage
- Carefully follow manufacturer's safety instructions on labels of chemical containers and material safety data sheets.



Caution

Feed water must be disconnected and the system must be depressurized before performing the cell cleaning procedure.

General Cleaning Instructions

Disconnect electrical service to the unit. Wipe exterior surfaces with lightly dampened cloth containing mild soap solution.

Cell Cleaning

1. Disconnect inlet water and depressurize. See "System Depressurization" section.
2. Turn unit OFF. Open right side door and remove top pc board cover.
3. Disconnect cable from J2 and snip plastic cable tie located on the inside top of unit.
4. Unscrew and remove the cell by turning the gray bushing counterclockwise. Carefully remove O-ring to clean the cell.
5. Wash the cell in a mild detergent solution or a 10% inorganic acid solution (follow acid manufacturer's recommended handling procedure). This may be done in an ultrasonic cleaner or with a soft brush.
6. Thoroughly rinse the cell in deionized or distilled water following the detergent or acid cleaning.
7. After cleaning, check the o-ring; replace if necessary.
8. Reinstall the cell into the cell well and hand tighten.
9. Reroute the cable and reconnect to J2 on circuit board.
10. Replace pc board cover and close door.
11. Reconnect water and turn unit ON.

**Note**

The sanitization timer will display after approximately six months, reminding you to sanitize the unit. Complete the system sanitization procedure and reset the sanitization timer according to the “Resetting the Sanitization Timer” section in this manual.

**Warning**

Disconnect from the power supply prior to maintenance and servicing.

**Warning**

Refer servicing to qualified personnel.

**Warning**

Avoid splashing disinfecting solution on clothing or skin.

Ensure all piping connections are tight to avoid leakage.

Ensure adequate ventilation.

Carefully follow manufacturer’s safety instructions and material safety data sheets.

This device to be used with water feeds only.

Sanitizing/cleaning agents must be used in compliance with instructions in this manual. Failure to comply with the above could result in explosion and personal injury.

System Sanitization

For UV/UF units, it is best to perform cleaning procedures early in the afternoon. After sanitization is complete (approximately 3.5 hours), place the unit into its normal (recirculation) mode to allow the system to recirculate overnight after cleaning.

1. From the (Idle) mode, press the UP or DOWN arrow until the display reads, “Sanitize Unit?”
2. Press ENTER.
3. Display will prompt, “Disconnected Inlet Water (NO).”
4. Disconnect feed water supply at the quick-disconnect inlet valve.
5. Press the UP or DOWN ARROW to select (YES) and then press ENTER. Display will now prompt, “Injected Sanitant (NO)?”
6. Remove solution syringe from wrapper.
7. Open front (left side) door. Remove luer cap on system injection port, located on the right side of the cartridge pack manifold by turning counterclockwise.
8. Remove the syringe luer cap and attach the syringe to the system luer fitting.
9. Slowly inject the solution into the system and remove syringe. Avoid injecting air into system.
10. Replace luer cap on injection port. Close door.
11. Press the UP or DOWN arrow to select (YES) and then press ENTER. Display will now read, “Reconnected Water (NO).”
12. Attach feed water supply removed in step 4.
13. Press the UP or DOWN arrow to select (YES) and press ENTER. Display will now read, “Unit Sanitizing.”

14. Unit may now be left unattended until the display reads, "Sanitization Complete: Press ENTER."
15. After sanitization is complete and the user presses "ENTER" the unit will return to the (Idle) mode.
16. Press the START/STOP switch to return the unit to the normal recirculating mode.

**Warning**

Depressurize system prior to removing the cartridge pack.

System Depressurization

This should be done prior to removing a cartridge pack or an ultrafilter. The unit needs to be plugged in and operating to correctly perform this depressurization step.

1. **IMPORTANT:** Disconnect inlet water from the left back of unit.
2. From the (Idle) display, press ENTER.
3. From "Air Purge?" press the UP arrow.
4. From "Depressurize" press the UP arrow to select (YES).
5. Press ENTER. Display will show "System Flushing" and the time remaining as the system depressurizes.

**Warning**

Depressurize system prior to removing the cartridge pack.

**Note**

A full cartridge pack may weigh about 20 lbs. (9 kg).

Cartridge Pack Replacement

The frequency with which you will need to clean your unit and replace your cartridge pack is dependent on your feed water's characteristics, your purity requirements and your usage. Sanitize your NANOpure Diamond TOC Life Science (UV/UF with TOC) and replace the cartridge pack when the product water purity drops below acceptable levels of resistivity, when organic levels become too high, or if a new 0.2 micron filter clogs rapidly after installation. The simple-to-use cleaning syringe is available from BarnsteadThermolyne (Catalog Number CMX25). This is used to effect a complete sanitization.

1. Remove the feed water line by depressing the stainless steel thumb pad.
2. Depressurize the system according to "System Depressurization" under the "Maintenance and Servicing" section.
3. Shut power OFF to the unit and disconnect the electrical service to the unit.
4. Open the front (left side) door.
5. Loosen the wing head screws on the cartridge manifold.
6. Remove the exhausted cartridge pack by lifting the unit manifold and pulling the cartridge pack out.
7. Install new cartridge pack, rinse and sanitize according to the instructions in the "Initial Operation" section.





Caution

Do not overtighten the 0.2 micron filter or use excessive force in seating it. The filter can be damaged by over-tightening or excessive force.



Note

If a newly installed 0.2 micron filter clogs rapidly after installation, the NANOpure Diamond TOC Life Science (UV/UF with TOC) may need to be sanitized to remove bacterial contaminants. See "System Sanitization."



0.2 Micron Filter



Warning

Depressurize the system before removing the ultrafilter.



Note

Some water will drain from the ultrafilter when it is disconnected. Make sure a sink, bucket or other water-proof container is available to place the ultrafilter in while it drains. Water will also drain from the tubing onto the bench or floor. To avoid slipping and falling, clean up all spills immediately.



Note

Do not install the ultrafilter until the NANOpure Diamond TOC Life Science (UV/UF with TOC) is ready for use. If the unit sits for longer than 96 hours without use, it will become contaminated and may affect product water quality.

0.2 Micron Filter Replacement

Replace the 0.2 micron final filter whenever any of the following conditions occur: every 30 days, the product water flow rate is reduced or bacteria break through. To replace the 0.2 micron filter assembly:

1. Remove the old 0.2 micron filter assembly by turning it counterclockwise until it is free from the dispense arm manifold.
2. Remove the new 0.2 micron filter assembly from the bag and insert it into the dispense arm manifold. Gently turn it clockwise until it is fully seated.
3. Remove the protective cap from the filter bell.
4. Rinse 1-2 liters of water through the filter to drain prior to using the product water.

Ultrafilter Replacement


The length of your ultrafilter's life will depend to some extent on conditions and use in your lab. When you are unable to completely sanitize your system (i.e. you can no longer get pyrogen-free water even after sanitizing the system), replace the ultrafilter as follows:

1. Disconnect feed water and depressurize the system according to "System Depressurization" under the "Maintenance and Servicing" section.
2. Turn the NANOpure Diamond TOC Life Science (UV/UF with TOC) OFF and disconnect it from the power supply.
3. Open the front (right side) door and right side door.
4. Separate the velcro strap from the ultrafilter.
5. Simultaneously push and hold the collet in using the tool supplied with the ultrafilter while pulling the top tubing out of the connector. Refer to photo on next page.
6. Unscrew tubing with luer fitting from the top side port of the ultrafilter.

MAINTENANCE AND SERVICING




Note
A small amount of water may drain from the ultrafilter when disconnected.

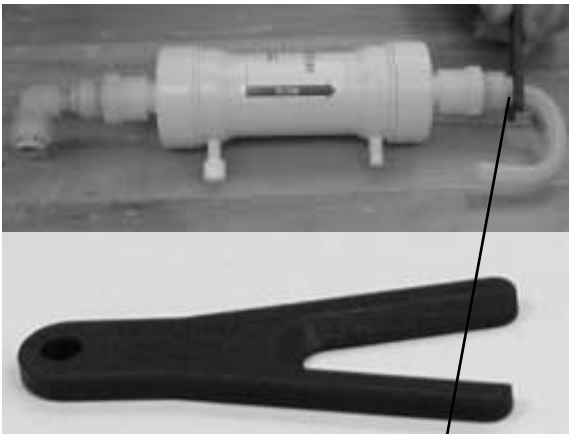


Tool used for ultrafilter removal

Collet

- 7. Using the tool supplied with the ultrafilter, simultaneously push up and hold the collet on the ultrafilter in while pulling the bottom fitting out of the connector. See photo.
- 8. Install new ultrafilter according to “Ultrafilter Installation” in the “Installation” section of this manual.
- 9. Sanitize the NANOpure Diamond TOC Life Science (UV/UF with TOC) according to the “System Sanitization” section of this manual.

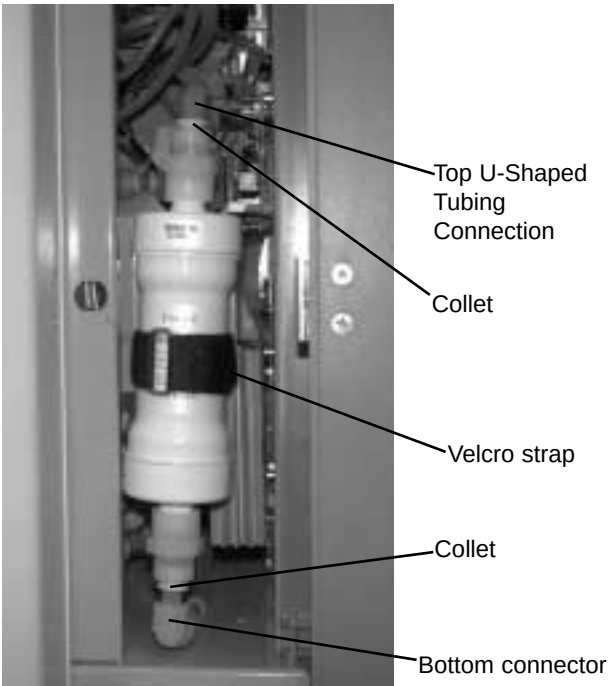

Ultrafilter Installation

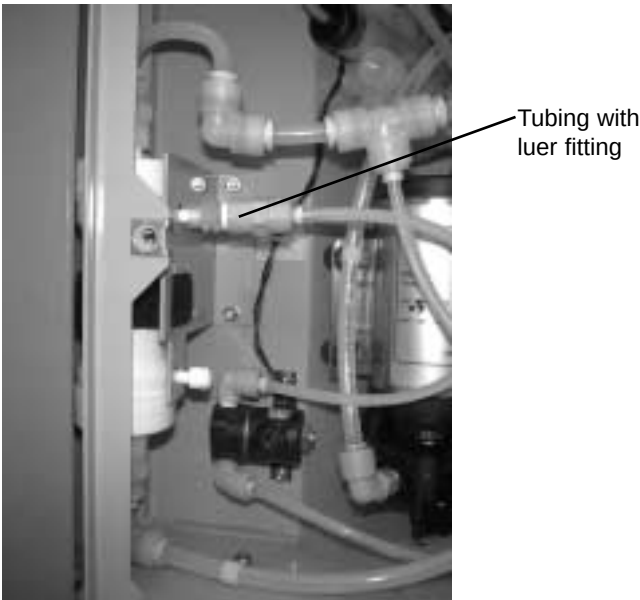

Ultrafilter Connections - Right side view



Warning

This unit is equipped with an ultraviolet lamp. Ultraviolet radiation is harmful to the eyes and skin. Do not attempt to observe the lamp directly.



Note

The NANOpure Diamond TOC Life Science (UV/UF with TOC) display will read "Check UV Lamp" when the lamp is burned out or disconnected.



Note

After approximately 5 months of operation, the display will periodically show a reminder to change the lamp.



Note

The UV lamp contains mercury. If broken or no longer needed, do not dispose of the UV lamp in the trash. Recycle or dispose of the UV lamp as hazardous waste.



UV Lamp Installation

UV Lamp Replacement

The ultraviolet lamp requires changing every 5-6 months. Lamp Life will vary according to the number of times the NANOpure Diamond TOC Life Science (UV/UF with TOC) is turned on and off (the START/STOP mode). This estimated lamp life is based on the NANOpure Diamond being operated in the normal recirculating mode during regular working hours and then placed in the Standby mode during off hours. Every fourth time the pump goes on in the Standby mode, the UV light will go on. If the NANOpure Diamond TOC Life Science (UV/UF with TOC) is cycled between the normal and Standby modes during the workday, this will result in a shorter lamp life. Therefore, it is recommended that the NANOpure Diamond (UV/UF with TOC) be left in the normal recirculating mode during regular working hours.

To replace the UV lamp:


1. Turn the NANOpure Diamond TOC Life Science (UV/UF with TOC) OFF and disconnect power cord from unit.
2. Access the left side of the unit by unlatching the screw securing the door with a flathead screwdriver. Unlatch and open the door. The UV lamp will be installed vertically.
3. Remove the black cap by sliding it off the UV chamber. Disconnect and remove old UV lamp and dispose of it in the proper manner (see Note).
4. Remove the UV lamp from its packaging. DO NOT TOUCH THE GLASS PORTION OF THE LAMP! It is recommended that lint-free gloves be worn when handling the lamp. The glass portion must be free of fingerprints, perspiration, etc. Even a light coating of perspiration will reduce the effectiveness of the lamp. If the glass portion of the lamp is touched, clean it with a damp, lint-free cloth: use isopropyl alcohol as required.
5. Insert the lamp halfway into the UV lamp chamber and plug UV cable in, then fully insert the

lamp into the chamber.

6. Make sure the UV lamp cable is tucked behind the chamber before closing and relatching the door.
7. Reconnect the power cord. Turn unit ON.
8. Reset UV Timer according to "Resetting the UV Timer."
9. Operate normally.

TOC Lamp Replacement

1. Place unit in the "Idle" mode.
2. Turn power off. Disconnect the power cord from the power entry module.
3. Open right side door and locate the TOC board at the bottom. The lamp chamber is protruding from the board on the right side.
4. Locate the clips on both sides of lamp chamber. While gently pulling out chamber, wedge a flat head screwdriver between chamber and clips to release chamber.
5. When installing the new lamp chamber, be sure to carefully insert pins of chamber into left side row of sockets on TOC board. See *Chamber Replacement* photo.
6. Finish installing the new lamp chamber by securing with clips. Clips may need to be pushed apart so the chamber will snap back into place.
7. Reconnect power cord and reset the TOC lamp timer. (See "Resetting the TOC Lamp Timer.")

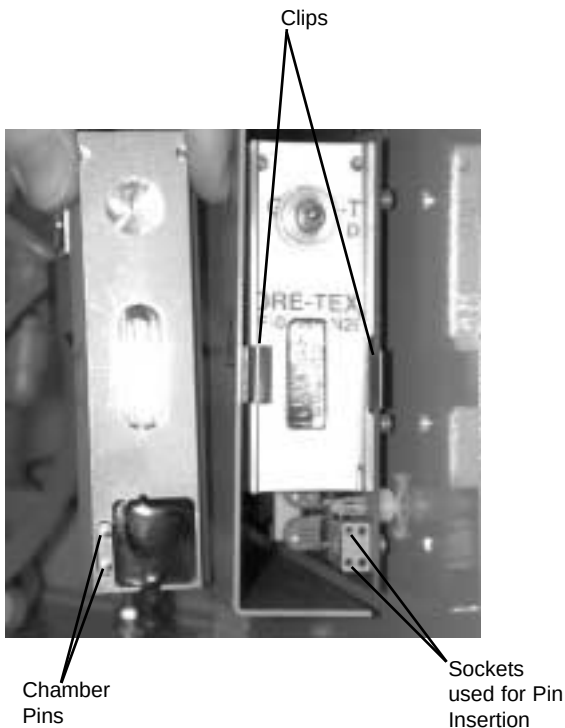

Chamber Replacement

**Warning**

Replace fuses with those of the same type and rating.

Fuse Replacement

1. Turn off the NANOpure Diamond TOC Life Science (UV/UF with TOC) and disconnect it from the power supply by removing the power cord directly below the main power switch on right side of unit. Locate the fuse drawer above the power cord receptacle.
2. Pull out the fuse drawer located in the power entry module.
3. Remove old fuses and replace with fuses of the same type and rating. (See "Replacement Parts.")
4. Replace fuse drawer.
5. Reattach the power cord.
6. Operate normally.

Shutdown

If the NANOpure Diamond TOC Life Science (UV/UF with TOC) is to be shut down for an extended period of time, the unit should be completely drained and the cartridge pack removed to prevent the growth of bacteria. If the system has remained inactive and full of water for more than 96 hours, the unit should be drained, sanitized and new cartridge pack installed prior to use.

Troubleshooting

Problem	Possible Causes	Solutions
NANOpure Diamond completely inactive (pump and display not operating).	No electrical power to NANOpure Diamond.	Ensure NANOpure Diamond power cord is connected to a live power source and completely plugged into an electrical outlet as well as power entry module on the unit. Make sure on/off switch is in the "i" (on) position.
	Main fuse(s) blown.	Replace the main fuse(s) as indicated in "Fuse Replacement."
	Power supply fuse blown.	Contact Customer Service concerning power supply fuse replacement.
Reduced or no product flow from the dispense arm or remote dispenser (if attached).	Red caps in cartridge pack.	Remove red caps from cartridge pack.
	0.2 micron filter clogged.	Replace final filter as indicated in "0.2 Micron Filter Replacement."
	UF module clogged.	Replace ultrafilter.
	Feed water not connected.	Reconnect feed water.
	Failed solenoid.	Replace BK1190X2.
Display reads, "Check Inlet" and pump does not run.	Pump protector (in reservoir), feed water line pressure switch or jumper plug not connected to pump interlock.	Connect the pump protector, jumper plug or pressure switch cord to the receptacle located on the upper left hand section of the Diamond.
	Insufficient feed water pressure or or volume to satisfy the NANOpure Diamond.	Open feed water line shut off, connect feed water line to Diamond and ensure sufficient pressure or sufficient water level in the feed water reservoir.
Display reads, "...MΩ-cm"	Air in cell.	If you just started your Diamond, wait 1 minute to see if it corrects itself.
	Bad cell connection.	Check connection of the cell as it attaches to the circuit board.
Recirculated water will not rinse up to desired purity level.	Exhausted cartridge pack.	Replace cartridge pack as indicated in "Cartridge Pack Replacement."

Problem	Possible Causes	Solutions
Short cartridge pack life.	Cartridge pack is beyond expiration date.	Check the expiration date. Cartridge packs begin to lose capacity after being stored for 2 years from the date of manufacture. Replace the cartridge pack with an unexpired one.
	Poor quality feed water.	If Barnstead ROpure is the feed water source, check that the membrane is functioning properly.
		If a Barnstead Still is the feed water source, ensure that the distillate temperature to the Diamond does not exceed 40°C (104°F).
		If feed water is from a central water purification system, verify water quality and proper functioning of the system. See “General Specifications.”
The display periodically reads, “Check UV Bulb.”	UV lamp not connected properly.	Check to ensure lamp is connected properly.
	UV lamp burned out.	Replace the UV lamp as described in “UV Lamp Replacement.”
Purity display continuously shows “OVERRANGE.”	Air in system.	Dispense more water through the NANOpure Diamond to purge air from the unit.
	Bad cell connection.	Check the connection of the cell as it attaches to the meter circuit board mounted on the inside. of the right door.
	Incorrect N.I.S.T. calibration offset.	If you have a N.I.S.T. module, perform a N.I.S.T. calibration.
		If you do not have a N.I.S.T. module, perform the calibration function WITHOUT the calibration module, as if you had it. (Unplug cell connection from J2 of PC1190X2 until finished.) Answer “Yes” to the prompt. The display will read “Calibration (failed),” but this will clear any incorrect offset to zero.
Display reads, “Overpressure Error.”	0.2 micron filter clogged.	Replace final filter as indicated in “0.2 Micron Filter Replacement.”
	Solenoid failure.	Replace BK1190X2

TROUBLESHOOTING

Problem	Possible Causes	Solutions
Display reads "Meter Board Problem."	PC board communication error.	Turn the unit OFF (I/O Switch) and restart.
	System electronics failure.	Replace the resistivity meter board mounted on the inside of the right door.
Water is leaking between the manifold and the cartridge pack.	Wing head screws are not tight enough.	Tighten wing head screws on unit manifold.
	New cartridge pack: o-ring missing from post (each of the 4 posts must contain an o-ring).	Disconnect inlet water and depressurize system. Next, disconnect the manifold and verify/replace missing o-ring. Reinstall cartridge pack as indicated in "Cartridge Pack Replacement."
	Debris on o-ring or bottom of unit manifold.	Inspect o-rings and bottom of unit manifold. Remove any objects with a lint-free cloth.
Upon system startup, the message "TOC Comm. Error: (0A)" is displayed.	PC board communication	Turn the Diamond OFF (I/O Switch) and restart the unit.
No TOC reading displayed.	TOC analyzer not activated.	Set "Use TOC Analyzer" to YES.
The display periodically reads "Replace TOC Bulb," but TOC timer has not reached zero.	TOC lamp is burned out or intensity is insufficient.	Replace TOC lamp. (See "TOC Lamp Replacement")
Display reads, "<TOC Error>"		
001	TOC memory error.	Cycle power, replace TOC if it continues.
002	TOC communications error.	Cycle power.
005, 007	TOC sensor error.	Cycle power, replace TOC if it continues.
006	Shorted TOC resistivity sensor.	Cycle power and place TOC analyzer in hour long clean mode (under Diagnostics); if still shorted, replace TOC.
All other error codes.		Contact technical customer service.

Replacement Parts

Recommended Spares

Consumables

Consumable parts are those required to support the day-to-day operation of this equipment.

Barnstead|Thermolyne establishes two types of consumables; those items that must periodically be replaced to maintain performance (filters, resin cartridges, etc.) and other items of limited life (indicator lights, fuses, etc.) that you can expect to replace on a more or less random basis. Where practical, Barnstead|Thermolyne recommends the frequency of replacement, or provides information on life expectancy from which you may calculate a replacement interval compatible with your usage pattern. The replacement of consumable parts is discussed in the Maintenance and Servicing section to assist you in accomplishing your own service.

Consumables may be ordered separately and in some cases, as an expendables kit. Check with your Barnstead|Thermolyne representative for additional information on the expendables kit.

<u>Description</u>	<u>Catalog No.</u>	<u>Recommended Quantity</u>
Diamond Kit Organic Free R/O & Distilled Feed	D50280	1
Diamond Kit Organic Free Deionized Feed	D50281	1
Final Filter	D3750	2
Main Fuse, 120 Volt - Slow Blow, 1.6 Amp, 250V	FZX47	2
Main Fuse, 230 Volt - Type T Time Lag, 0.63 Amp, 250V	FZX54	2
Teflon [®] Tape, Roll	06078	1
UV Lamp	LMX13	2
Transformer Fuse, Type F Fast-Acting, 3.15 Amp, 250 V	AZ9027	1
Sanitization Solution	CMX25	2
Ultrafilter	FL1192X1	1
TOC Analyzer Lamp	LMX26	1

General Maintenance Parts

General maintenance parts are defined as laboratory level repair parts which do not require great expertise or special tools for installation. **Barnstead|Thermolyne** recommends that you stock the general maintenance parts as an aid to ensuring the continued operation of this equipment.

<u>Description</u>	<u>Catalog No.</u>	<u>Recommended Quantity</u>
Quick Disconnect Inlet	CUX11	1
Quick Disconnect Insert 3/8" tube	CUX9	1
Cartridge Pack O-Rings	GSX68	4

[®]Teflon is a registered trademark of DuPont.

REPLACEMENT PARTS

Safety Stock


For critical applications where performance with minimum downtime is required, **Barnstead|Thermolyne** recommends that you maintain a local stock of those parts listed in the General Maintenance Parts and Safety Stock sections.

<u>Description</u>	<u>Part No.</u>	<u>Recommended Quantity</u>
Display	SW1190X2	1
External Display cable (Remote Display Use)	WHX20	1
Main Control Board (PC1)	**See page 54.	
Recirculation Pump & Motor	PU1190X1	1
Pressure Regulator	PM1190X1	1
Resistivity Cell	E896X4	1
Solid State Ballast UV	SC1191X1	1
Quartz Sleeve, UV Lamp	TU733X1	1
Power Supply PC Board	TNX116	1
Calibration PC Board	PC1190X2	1
Resistivity Meter PC Board	PC896X2	1
Check Valve	PMX214	1
Wire Harness, UV Bulb	WH1191X1	1
Flow Sensor	ME1190X1	1
Flush Solenoid	RY1190X3	1
Cartridge Pack Manifold	BK1190X3	1
Dispenser Manifold W/Solenoids	BK1190X2	1
Dispense Switch	SW1190X1	1
Cell O-ring	GSX29	1
TOC Analyzer	PCX54	1

Optional Accessory Ordering Information

<u>Description</u>	<u>Catalog No.</u>
N.I.S.T. Calibration Module	E896X5
Low-Pressure Switch	D2706
Float Switch	D8964
Remote Dispenser	D11981
Serial Cable for Communication with Computer	WHX18
Printer	AY1137X1
Printer Paper	AY669X4

Flow Chart

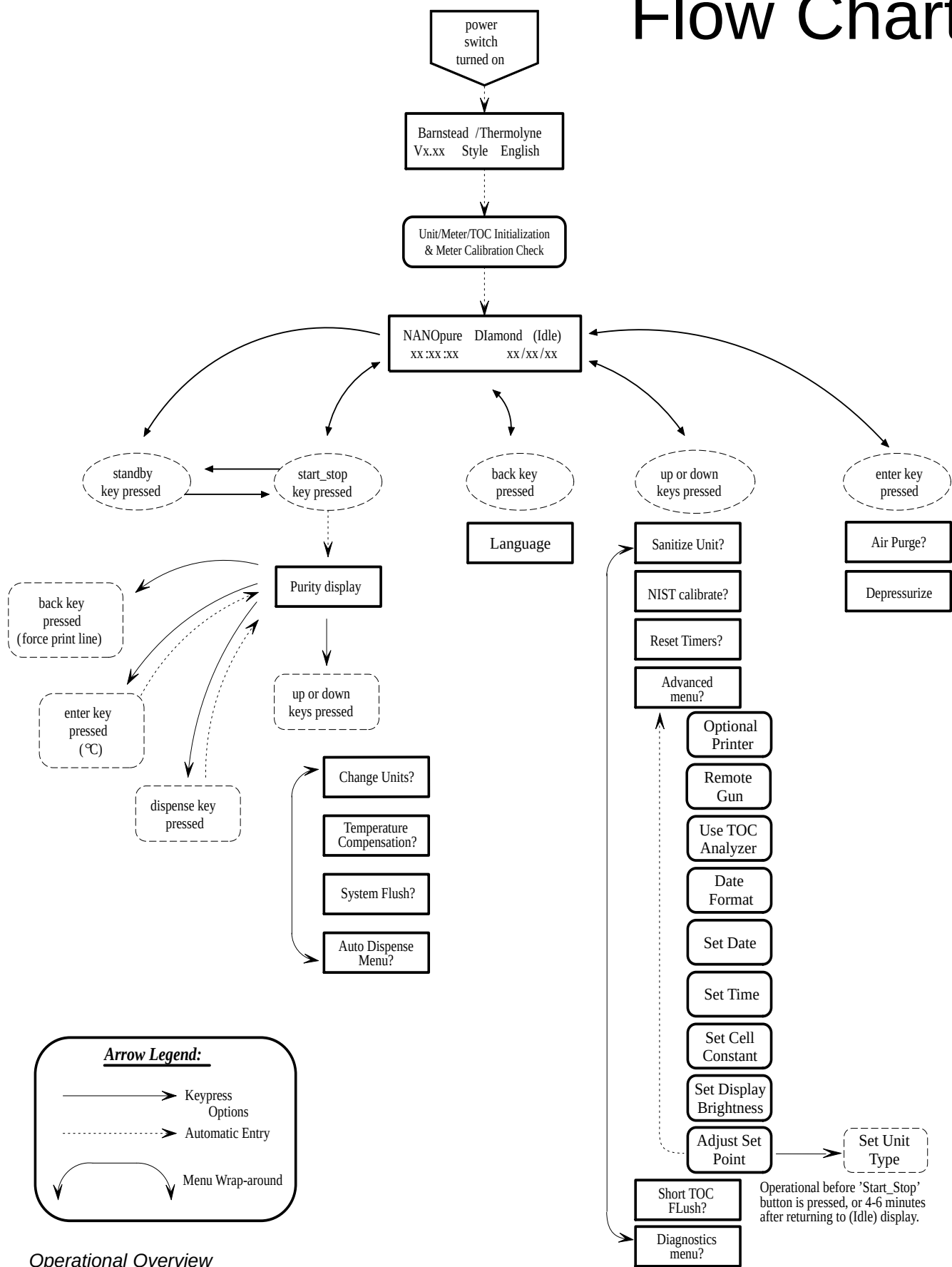

Operational Overview

Wiring Diagram

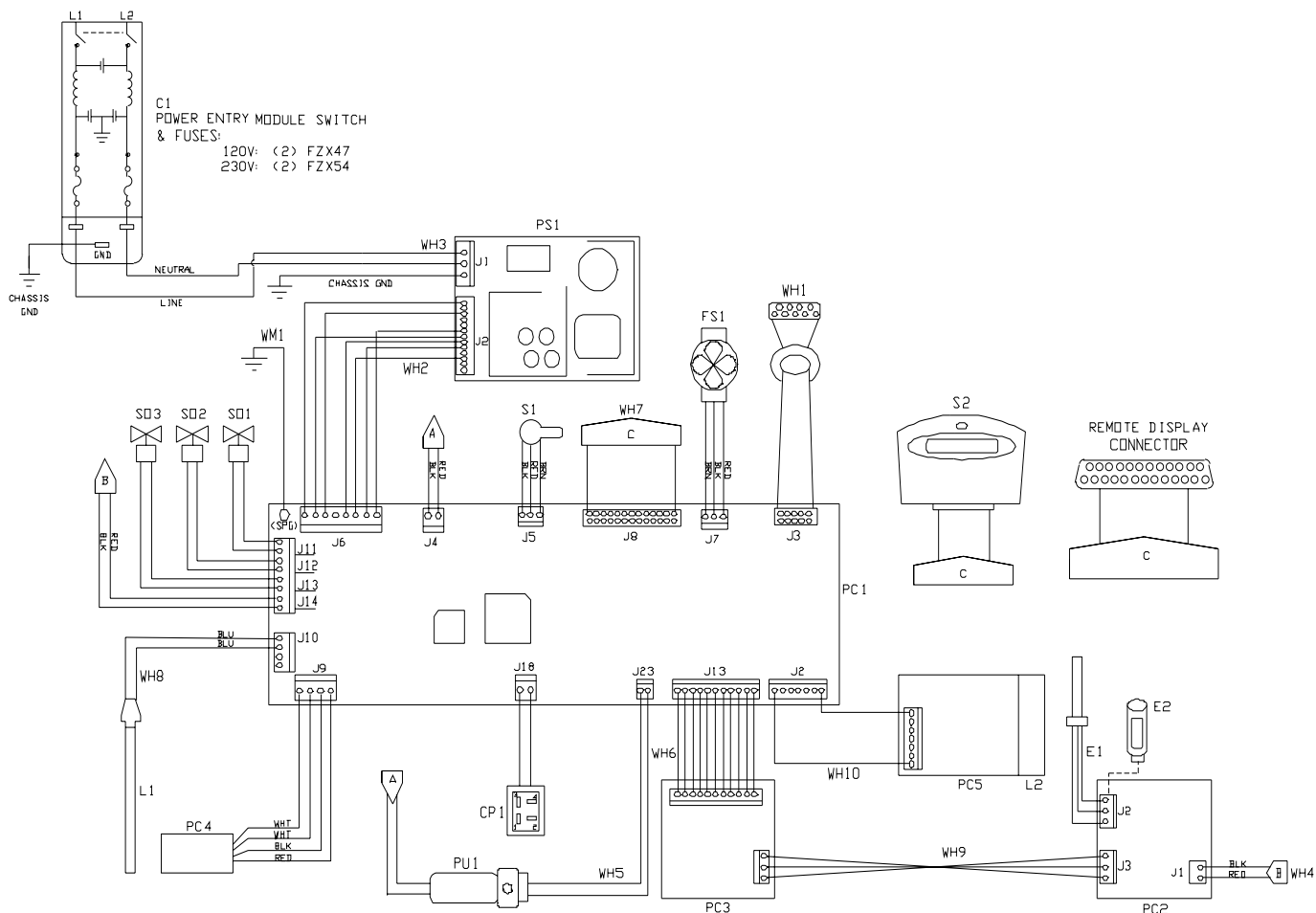


DIAGRAM COMPONENT LIST

REF. NO.	DESCRIPTION	MODEL NO. AND OUR PART NO. (S)											
		D11901	D11911	D11921	D11931	D11941	D11951	D11961	D11971	D12411	D12421	D12431	D12441
C1	POWER ENTRY MODULE	CEX183	CEX183	CEX183	CEX183	CEX183	CEX183	CEX183	CEX183	CEX183	CEX183	CEX183	CEX183
CP1	PUMP INTERLOCK	CEB96X2	CEB96X2	CEB96X2	CEB96X2	CEB96X2	CEB96X2	CEB96X2	CEB96X2	CEB96X2	CEB96X2	CEB96X2	CEB96X2
E1	CONDUCTIVITY CELL	E9B6X4	E9B6X4	E9B6X4	E9B6X4	E9B6X4	E9B6X4	E9B6X4	E9B6X4	E9B6X4	E9B6X4	E9B6X4	E9B6X4
E2	NIST MODULE	E9B6X5	E9B6X5	E9B6X5	E9B6X5	E9B6X5	E9B6X5	E9B6X5	E9B6X5	E9B6X5	E9B6X5	E9B6X5	E9B6X5
FS1	FLOW SENSOR	ME1190X1	ME1190X1	ME1190X1	ME1190X1	ME1190X1	ME1190X1	ME1190X1	ME1190X1	ME1190X1	ME1190X1	ME1190X1	ME1190X1
L1	UV LAMP	—	LMX13	—	—	LMX13	—	—	LMX13	—	LMX13	—	LMX13
L2	TOC LAMP ASSEMBLY	—	—	—	—	LMX26	LMX26	LMX26	LMX26	LMX26	LMX26	LMX26	LMX26
PC1	MAIN CONTROL BOARD	PC1190X4	PC1191X1	PC1192X1	PC1193X1	PC1194X2	PC1195X1	PC1196X1	PC1197X1	PC1241X2	PC1242X1	PC1243X1	PC1244X1
PC2	CALIBRATION BOARD	PC1190X2	PC1190X2	PC1190X2	PC1190X2	PC1190X2	PC1190X2	PC1190X2	PC1190X2	PC1190X2	PC1190X2	PC1190X2	PC1190X2
PC3	METER BOARD	PCB96X2	PCB96X2	PCB96X2	PCB96X2	PCB96X2	PCB96X2	PCB96X2	PCB96X2	PCB96X2	PCB96X2	PCB96X2	PCB96X2
PC4	BALLAST	—	SC1191X1	—	—	—	SC1191X1	—	SC1191X1	—	SC1191X1	—	SC1191X1
PC5	TOC ANALYZER	—	—	—	—	PCX54	PCX54	PCX54	PCX54	PCX54	PCX54	PCX54	PCX54
PU1	PUMP	PU1190X1	PU1190X1	PU1190X1	PU1190X1	PU1190X1	PU1190X1	PU1190X1	PU1190X1	PU1190X1	PU1190X1	PU1190X1	PU1190X1
PS1	POWER SUPPLY	TNX116	TNX116	TNX116	TNX116	TNX116	TNX116	TNX116	TNX116	TNX116	TNX116	TNX116	TNX116
S1	DISPENSE SWITCH	SW1190X1	SW1190X1	SW1190X1	SW1190X1	SW1190X1	SW1190X1	SW1190X1	SW1190X1	SW1190X1	SW1190X1	SW1190X1	SW1190X1
S2	MEMBRANE SWITCH	SW1190X2	SW1190X2	SW1190X2	SW1190X2	SW1190X2	SW1190X2	SW1190X2	SW1190X2	SW1190X2	SW1190X2	SW1190X2	SW1190X2
SD1	DIVERTER SOLENOID	RY1190X2	RY1190X2	RY1190X2	RY1190X2	RY1190X2	RY1190X2	RY1190X2	RY1190X2	RY1190X2	RY1190X2	RY1190X2	RY1190X2
SD2	DRAWOFF SOLENOID	RY1190X1	RY1190X1	RY1190X1	RY1190X1	RY1190X1	RY1190X1	RY1190X1	RY1190X1	RY1190X1	RY1190X1	RY1190X1	RY1190X1
SD3	FLUSH SOLENOID	RY1190X3	RY1190X3	RY1190X3	RY1190X3	RY1190X3	RY1190X3	RY1190X3	RY1190X3	RY1190X3	RY1190X3	RY1190X3	RY1190X3
WH1	PRINTER HARNESS	WH1190X1	WH1190X1	WH1190X1	WH1190X1	WH1190X1	WH1190X1	WH1190X1	WH1190X1	WH1190X1	WH1190X1	WH1190X1	WH1190X1
WH2	DC POWER HARNESS	WH1190X2	WH1190X2	WH1190X2	WH1190X2	WH1190X2	WH1190X2	WH1190X2	WH1190X2	WH1190X2	WH1190X2	WH1190X2	WH1190X2
WH3	AC POWER HARNESS	WH1190X3	WH1190X3	WH1190X3	WH1190X3	WH1190X3	WH1190X3	WH1190X3	WH1190X3	WH1190X3	WH1190X3	WH1190X3	WH1190X3
WH4	CARD BOARD HARNESS	WH1190X4	WH1190X4	WH1190X4	WH1190X4	WH1190X4	WH1190X4	WH1190X4	WH1190X4	WH1190X4	WH1190X4	WH1190X4	WH1190X4
WH5	PRESSURE SW. HARNESS	WH1190X5	WH1190X5	WH1190X5	WH1190X5	WH1190X5	WH1190X5	WH1190X5	WH1190X5	WH1190X5	WH1190X5	WH1190X5	WH1190X5
WH6	METER BOARD HARNESS	WH1190X6	WH1190X6	WH1190X6	WH1190X6	WH1190X6	WH1190X6	WH1190X6	WH1190X6	WH1190X6	WH1190X6	WH1190X6	WH1190X6
WH7	DISPLAY HARNESS	WXH19	WXH19	WXH19	WXH19	WXH19	WXH19	WXH19	WXH19	WXH19	WXH19	WXH19	WXH19
WH8	LAMP HARNESS	—	WH1191X1	—	WH1191X1	—	WH1191X1	—	WH1191X1	—	WH1191X1	—	WH1191X1
WH9	CH / METER HARNESS	WM896X2A	WM896X2A	WM896X2A	WM896X2A	WM896X2A	WM896X2A	WM896X2A	WM896X2A	WM896X2A	WM896X2A	WM896X2A	WM896X2A
WH10	TOC HARNESS	—	—	—	—	WH1194X1	WH1194X1	WH1194X1	WH1194X1	WH1194X1	WH1194X1	WH1194X1	WH1194X1
WM1	GROUND WIRE	WM1190X2	WM1190X2	WM1190X2	WM1190X2	WM1190X2	WM1190X2	WM1190X2	WM1190X2	WM1190X2	WM1190X2	WM1190X2	WM1190X2

Ordering Procedures

Please refer to the Specification Plate located rear right-hand side for the complete model number, serial number, and series number when requesting service, replacement parts or in any correspondence concerning this unit.

All parts listed herein may be ordered from the Barnstead|Thermolyne dealer from whom you purchased this unit or can be obtained promptly from the factory. When service or replacement parts are needed we ask that you check first with your dealer. If the dealer cannot handle your request, then contact our Customer Service Department at 319-556-2241 or 800-553-0039.

Prior to returning any materials to Barnstead|Thermolyne Corp., please contact our Customer Service Department for a "Return Goods Authorization" number (RGA). Material Returned without an RGA number will be returned.

Two Year Limited Warranty

BARNSTEAD|THERMOLYNE CORPORATION ("BARNSTEAD") warrants that if a product manufactured by Barnstead shall be free of defects in materials and workmanship for two (2) year from the first to occur of (i) the date the product is sold by BARNSTEAD or (ii) the date the product is purchased by the original retail customer (the "Commencement Date"). Except as expressly stated above, BARNSTEAD MAKES NO OTHER WARRANTY, EXPRESSED OR IMPLIED, WITH RESPECT TO THE PRODUCTS AND EXPRESSLY DISCLAIMS ANY AND ALL WARRANTIES, INCLUDING BUT NOT LIMITED TO, WARRANTIES OF DESIGN, MERCHANT ABILITY AND FITNESS FOR A PARTICULAR PURPOSE.

An authorized representative of BARNSTEAD must perform all warranty inspections. In the event of a defect covered by BARNSTEAD's warranty, BARNSTEAD shall, as its sole obligation and exclusive remedy, provide free replacement parts to remedy the defective product. In addition, for products sold by BARNSTEAD within the continental United States or Canada, BARNSTEAD shall provide provide free labor to repair the products with the replacement parts, but only for a period of ninety (90) days from the Commencement Date.

BARNSTEAD's warranty provided hereunder shall be null and void and without further force or effect if there is any (i) repair made to the product by a party other than BARNSTEAD or its duly authorized service representative, (ii) misuse (including use inconsistent with written operating instructions for the product), mishandling, contamination, overheating, modification or alteration of the product by any customer or third party or (iii) use of replacement parts that are obtained from a party who is not an authorized dealer of BARNSTEAD.

Heating elements, because of their susceptibility to overheating and contamination, must be returned to the BARNSTEAD factory and if, upon inspection, it is concluded that failure is due to factors other than excessive high temperature or contamination, BARNSTEAD will provide warranty replacement. As a condition to the return of any product, or any constituent part thereof, to BARNSTEAD's factory, it shall be sent prepaid and a prior written authorization from BARNSTEAD assigning a Return Goods Number to the product or part shall be obtained.

IN NO EVENT SHALL BARNSTEAD BE LIABLE TO ANY PARTY FOR ANY DIRECT, INDIRECT, SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES, OR FOR ANY DAMAGES RESULTING FROM LOSS OF USE OR PROFITS, ANTICIPATED OR OTHERWISE, ARISING OUT OF OR IN CONNECTION WITH THE SALE, USE OR PERFORMANCE OF ANY PRODUCTS, WHETHER SUCH CLAIM IS BASED ON CONTRACT, TORT (INCLUDING NEGLIGENCE), ANY THEORY OF STRICT LIABILITY OR REGULATORY ACTION.

The name of the authorized Barnstead|Thermolyne dealer nearest you may be obtained by calling 1-800-446-6060 (319-556-2241) or writing to:

Barnstead|Thermolyne
P.O. Box 797
2555 Kerper Boulevard
Dubuque, IA 52004-0797
USA
FAX: (319) 589-0516
E-MAIL ADDRESS: mkt@barnstead.com